


AMT Bridle™

Family of Nasal Tube Retaining Systems


PEDIATRIC + ADULT

Placement Reference Guide


This is not a substitute for the DFU/IFU (Directions for Use/Instructions for Use). ⓘ ⓘ

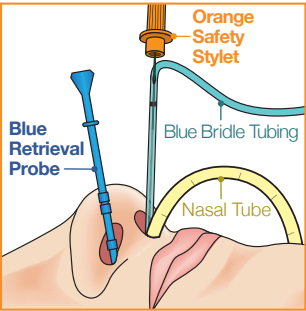


ADVANCE
Advance the blue retrieval probe into the nostril without the nasal tube. On the opposite side, insert the blue tubing (with orange safety stylet inside) until the magnets connect.

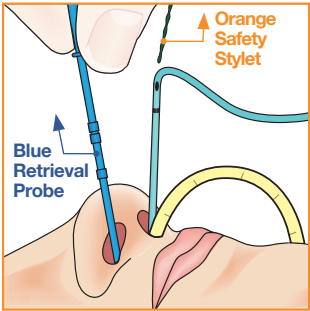


REMOVE
First, remove the orange safety stylet. Second, remove the blue retrieval probe, drawing the blue bridge tubing around the vomer bone and then out the opposite nare.

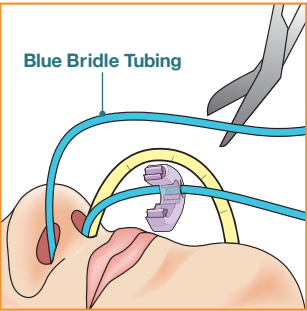

CUT
Cut and dispose of the blue bridge tube section containing the magnet. Slide the clip up the blue bridge tubing into position, and ensure the clip opens medially.

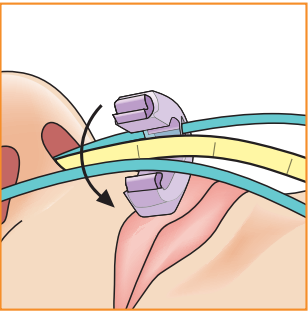


SECURE
Secure the nasal tube in the clip's groove. Place the blue bridge tubing within the clip (tubing position varies between clip styles).

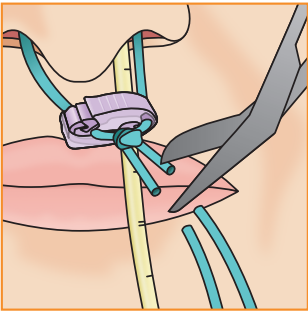


CLOSE
Close the clip 1 cm (1/2") below the nose. Tie 2-3 simple knots and cut off any excess.

Vomer Bone Securement

Patient Safety + Comfort are a top priority for AMT and the vomer bone is a great securement location for nasal tubes. Once an AMT Bridle™ Family device is properly placed, no tape or sutures are required to keep the nasal tube secured.

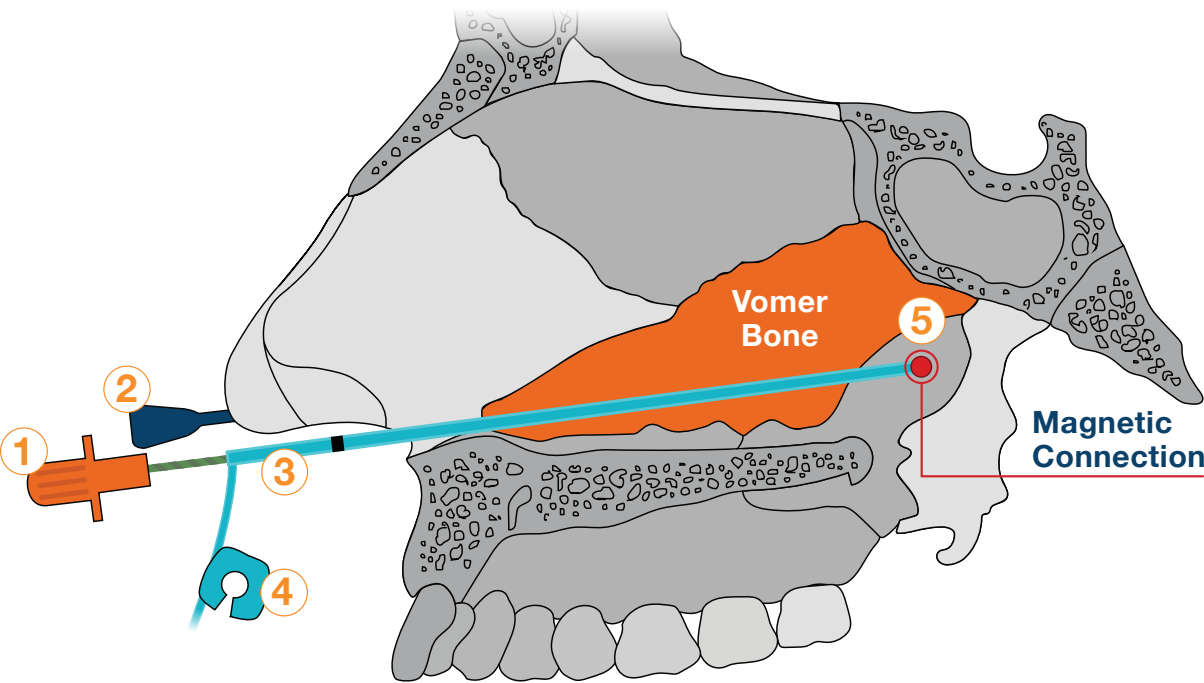


- 1 Safety Stylet
- 2 Retrieval Probe
- 3 Bridle Tubing
- 4 Pre-Attached Clip
- 5 Rare Earth Magnets
- 6 Clip Opening Tool

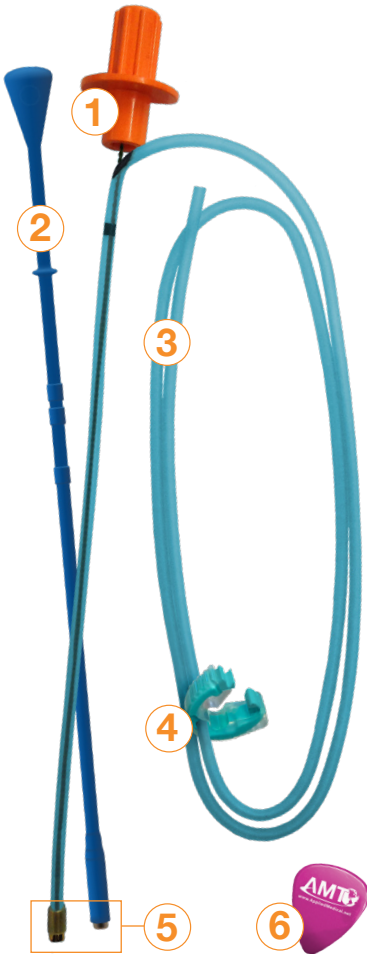

Accommodates Salem Sump®

The Bridle Pro® 14F, 16F and 18F clips, with specially designed openings, can accommodate large bore nasal tubes, including the Salem Sump®.

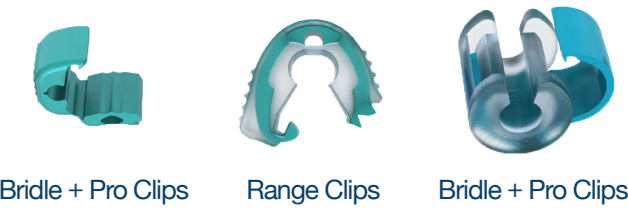

Bridle Family | Size Availability*

	Pediatric Bridle Pro® Range¹	Pediatric Bridle Pro®¹	Bridle Pro® Range¹	Bridle Pro®¹	AMT Bridle™²
5F	4-420506P				
6F					
8F	4-420810P	4-4208P	4-420810	4-4208	4-4108
10F				4-4210	4-4110
12F			4-421214	4-4212	4-4112
14F				4-4214	4-4114
16F			4-421618	4-4216	4-4116
18F				4-4218	4-4118

The Pro Range Clip's flexible overmolded section allows it to encompass multiple sizes into a single clip. ¹The Bridle Pro® features the exclusive blue Bridle Tubing. | ²The AMT Bridle™ features white umbilical tape. | *The Bridle Family comes packaged as 5/Box.

Questions? Call Today:

Jälleenmyyjä:
MEDI PLAST
Medioplast Fenno Oy
PL 153, 01531 Vantaa
Puh: 09 276 360
asiakaspalvelu@mediplast.com
www.mediplast.com