



# danumed® Extension Set for Gastrostomy Button

*Extension set for delivery of enteral nutrition, liquid medication and decompression.*

*Innovative features optimise the flow rate and simplify the handling for the user.*

*Cost-efficient solution with usage for up to 14 days.*



## COMFORT




### Innovative Rotating Connector ENSwivel

- Allows for the easy connection and release of ENFit® connections.
- Even very tight connections can be easily released.
- Only the rotating Swivel is turned and not the entire connector. This feature **avoids any twisting of the enteral tube** which is both painful for the patient and inconvenient for the user.
- All danumed ENFit® connectors offer the advantages of the ENSwivel®: **ENFit®**, **Y-ENFit®** and **Bolus-ENFit®**.

### High-End Button Connector

- The locking mechanism with optimised design ensures **easy, safe and tight connection** to the button.

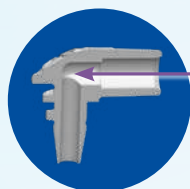

The rounded kink in the inner lumen optimises the flow rate, independent of the tube lumen.



Ergonomic form of the button connector with ribbed profile for a better grip.

## SAFETY

### Maximum safety for the patient

- The ENFit® connector according to DIN EN ISO 80369-3 **avoids tubing misconnections** and enhances patient safety significantly.
- All ENFit® connectors from danumed are made of **TRITAN™**, a plastic material resistant to aggressive solutions/medicines.
- All extension sets with ENFit® connector are **free from latex, PVC, DEHP and BPA**.
- The biocompatible tube made of polyurethane is **kink resistant** and **reduces the hazard of tube occlusions** and the number of alarm signals of the enteral pump.



## danumed® Extension Set with angled Button Connector



REF	Description	Length (cm)	Feeding connector	Latex-free	PVC-free	DEHP-free	BPA-free
DA 1011 3011	Extension Set <b>sterile</b> , 30 cm, ENFit®, ENSwivel®	30 cm	ENFit®	•	•	•	•
DA 1010 3011	Extension Set <b>non-sterile</b> , 30 cm, ENFit®, ENSwivel®	30 cm	ENFit®	•	•	•	•
DA 1011 6011	Extension Set <b>sterile</b> , 60 cm, ENFit®, ENSwivel®	60 cm	ENFit®	•	•	•	•
DA 1010 6011	Extension Set <b>non-sterile</b> , 60 cm, ENFit®, ENSwivel®	60 cm	ENFit®	•	•	•	•
DA 1011 0521	Extension Set <b>sterile</b> , 5 cm, Y-ENFit®, ENSwivel®	5 cm	Y-ENFit®	•	•	•	•
DA 1011 3021	Extension Set <b>sterile</b> , 30 cm, Y-ENFit®, ENSwivel®	30 cm	Y-ENFit®	•	•	•	•
DA 1010 3021	Extension Set <b>non-sterile</b> , 30 cm, Y-ENFit®, ENSwivel®	30 cm	Y-ENFit®	•	•	•	•
DA 1011 6021	Extension Set <b>sterile</b> , 60 cm, Y-ENFit®, ENSwivel®	60 cm	Y-ENFit®	•	•	•	•
DA 1010 6021	Extension Set <b>non-sterile</b> , 60 cm, Y-ENFit®, ENSwivel®	60 cm	Y-ENFit®	•	•	•	•

Delivered with right-angled button connector and **standard bore tube**. Individually packaged. 10 pieces in a box, 10 boxes in a shipping box.



## danumed® Extension Set with straight Button Connector



REF	Description	Length (cm)	Feeding connector	Latex-free	PVC-free	DEHP-free	BPA-free
DA 1011 3032	Extension Set <b>sterile</b> , 30 cm, Bolus-ENFit®, ENSwivel®	30 cm	Bolus-ENFit®	•	•	•	•
DA 1010 3032	Extension Set <b>non-sterile</b> , 30 cm, Bolus-ENFit®, ENSwivel®	30 cm	Bolus-ENFit®	•	•	•	•
DA 1011 6032	Extension Set <b>sterile</b> , 60 cm, Bolus-ENFit®, ENSwivel®	60 cm	Bolus-ENFit®	•	•	•	•
DA 1010 6032	Extension Set <b>non-sterile</b> , 60 cm, Bolus-ENFit®, ENSwivel®	60 cm	Bolus-ENFit®	•	•	•	•

Delivered with straight button connector and **wide bore tube**. Individually packaged. 10 pieces in a box, 10 boxes in a shipping box. The Bolus Extension Sets are available on request with funnel connector instead of ENFit® connector.

## Related Products to danumed® Extension Sets:



danuButton®  
Gastrostomy Button



danumed®  
Enteral Syringes, Single Use  
Enteral Syringes, Single-Patient Use

**danumed Medizintechnik GmbH –**  
Individual solutions are our standard.

We are an innovative, German medical technology company with the highest quality standards, many years of expertise in enteral nutrition and clear customer focus. We also meet individual customer requests with flexibility and reliability. Personal, competent, collaborative.

danumed Medizintechnik GmbH • Stadthof 18 • 93059 Regensburg • Germany  
Tel + 49 (0) 941/ 890 588-80 • Fax + 49 (0) 941/ 890 588-81

info@danumed.com • [www.danumed.com](http://www.danumed.com)

Forhandler:  
**MEDI PLAST**  
Medioplast AS  
Bjørnstjerne Bjørnsonsgate 110  
3044 Drammen  
T 32 88 11 00  
[mediplast@mediplast.no](mailto:mediplast@mediplast.no)  
[www.mediplast.com](http://www.mediplast.com)

