

HARRY HOLMS

PRODUCT OVERVIEW

Harry Holms is the innovative specialist in stainless steel products. We offer an appreciated range of wire baskets, instrument trays, trolleys, holders and dispensers.



HARRY HOLMS™
STAINLESS STEEL PRODUCTS SWEDEN

STAINLESS STEEL TROLLEY

Harry Holms AB, in consultation with the majority of Swedish hospitals, has developed a stainless stacking trolley designed to be equipped with autoclave baskets, grid shelves with or without edges intended for containers or other goods.

The design patent guides provide easy handling and good fixation of the mentioned accessories.

Autoclave baskets and grid shelves for containers can be pulled out to half their depth from both sides of the trolley. The guides for autoclave baskets provide a slightly downward angle for easier loading and unloading of goods. The trolley can be loaded from both sides.

The trolley is designed for efficient cleaning without parts that create dirt pockets and consists only of stainless steel details.

The rounded top construction and the absence of protruding details provide smooth handling in connection with covering the carriage.

The four large (Ø125 mm) and lockable quality wheels create a prerequisite for a quiet and compliant cart. The wheels can withstand a temperature of up to 95°C.



STAINLESS STEEL TROLLEY

Art.no.	Description
883917C	1 trolley, 2 shelves 3 baskets



HARRY HOLMS™
STAINLESS STEEL PRODUCTS SWEDEN

INSTRUMENT TRAY, SPRI

Instrument tray intended to be used for cleaning and sterilizing of various surgical instruments and tools. The tray serves as a carrier/holder of instruments and tools during internal transport and during cleaning/sterilization in various dish disinfectors and autoclaves.

The product is a Medical device of risk class I (MDR EU 2017/74)

Material: Stainless steel EN 1.4301 C/C

Grid size: 8 mm

Handle the product carefully so that it is not subjected to shocks or impacts



INSTRUMENT TRAY, SPRI

Art.no.	UDI-DI	Description	Measure	Max load
705010	735012544003	SPRI III	250x170x70 mm	14 KG
705005	735012544002	SPRI II	340x250x70 mm	12 KG
705000	735012544001	SPRI I	450x340x70 mm	10 KG



HARRY HOLMS™
STAINLESS STEEL PRODUCTS SWEDEN

LID TO INSTRUMENT TRAY, SPRI


Lid intended to be used, among other things, for cleaning and sterilizing of various surgical instruments and tools. The lid functions as a closure of grids for instruments and tools during internal transport and during cleaning/sterilisation in various dish disinfectors and autoclaves.

The product is a Medical device of risk class I (MDR EU 2017/74)

Material: Stainless steel EN 1.4301 C/C

Grid size: 8 mm

Handle the product carefully so that it is not subjected to shocks or impacts

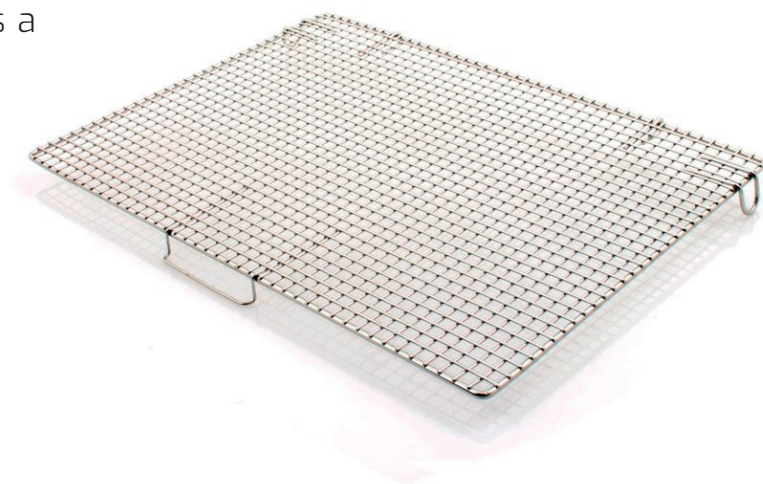

LID TO INSTRUMENT TRAY, SPRI

Art.no.	UDI-DI	Description	Measure
705210	735012544006	SPRI III	250x170 mm
705205	735012544005	SPRI II	340x250 mm
705200	735012544004	SPRI I	450x340 mm



HARRY HOLMS™
STAINLESS STEEL PRODUCTS SWEDEN

INNER LID TO INSTRUMENT TRAY, SPRI

Inner lid intended to be used, among other things, for cleaning and sterilizing of various surgical instruments and tools. The lid functions as a carrier/holder of instruments and tools during internal transport and during cleaning/sterilization in different dish disinfectors and autoclaves.



The product is a Medical device of risk class I (MDR EU 2017/74)

Material: Stainless steel EN 1.4301 C/C

Grid size: 8 mm

Handle the product carefully so that it is not subjected to shocks or impacts

INNER LID TO INSTRUMENT TRAY, SPRI			
Art.no.	UDI-DI	Description	Measure
705225	735012544009	SPRI III	237x157 mm
705220	735012544008	SPRI II	327x237 mm
705215	735012544007	SPRI I	437x327 mm



HARRY HOLMS™
STAINLESS STEEL PRODUCTS SWEDEN

INSTRUMENT TRAY, FINE MESH


Instrument tray intended to be used for cleaning and sterilizing of various surgical instruments and tools. The tray serves as a carrier/holder of instruments and tools during internal transport and during cleaning/sterilization in various dish disinfectors and autoclaves.

The product is a Medical device of risk class I (MDR EU 2017/74)

Material: Stainless steel EN 1.4301 C/C

Grid size: 1.6 mm

Handle the product carefully so that it is not subjected to shocks or impacts

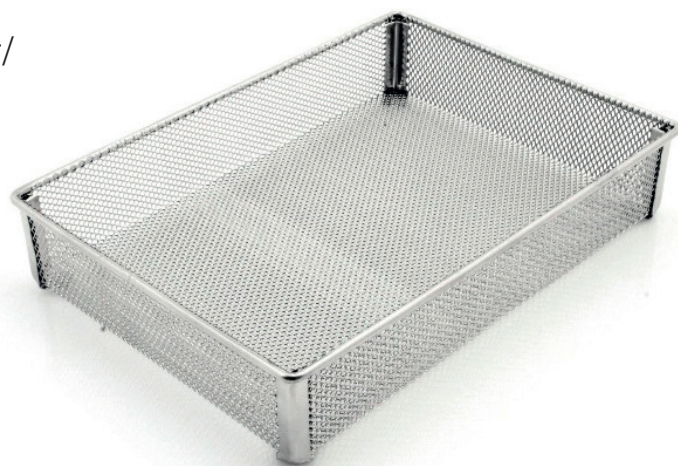

INSTRUMENT TRAY, FINE MESH

Art.no.	UDI-DI	Measure	Max load
705300	735012544010	420x310x40 mm	8 KG
705305	735012544011	310x210x40 mm	10 KG
705310	735012544012	210x150x40 mm	12 KG



HARRY HOLMS™
STAINLESS STEEL PRODUCTS SWEDEN

LID TO INSTRUMENT TRAY, FINE MESH

Lid intended to be used, among other things, for cleaning and sterilizing of various surgical instruments and tools. The lid functions as a closure of trays for instruments and tools during internal transport and during cleaning/sterilization in various dish disinfectors and autoclaves.

The product is a Medical device of risk class I (MDR EU 2017/74)

Material: Stainless steel EN 1.4301 C/C

Grid size: 1.6 mm

Handle the product carefully so that it is not subjected to shocks or impacts



LID TO INSTRUMENT TRAY, FINE MESH			
Art.no.	UDI-DI	Measure	Max load
705330	735012544015	210x150 mm	10 KG
705325	735012544014	310x210 mm	8 KG
705320	735012544013	420x310 mm	6 KG



HARRY HOLMS™
STAINLESS STEEL PRODUCTS SWEDEN

INSTRUMENT TUBE WITH LID

Washing basket intended to be used, among other things, for cleaning and sterilizing various surgical instruments and tools. The basket functions as a carrier/holder for instruments and tools during internal transport and during cleaning/sterilization in various dish disinfectors and autoclaves.

The product is a Medical device of risk class I (MDR EU 2017/74)

Material: Stainless steel EN 1.4301 C/C
Grid size: 1.6 mm



INSTRUMENT TUBE WITH LID

Art.no.	UDI-DI	Measure
705365	735012544031	Ø48x403 mm
705364	735012544024	Ø48x218 mm



HARRY HOLMS™
STAINLESS STEEL PRODUCTS SWEDEN

INSTRUMENT BASKET WITH LID

Washing basket intended to be used, among other things, for cleaning and sterilizing of various surgical instruments and tools.

The basket functions as a carrier/holder for instruments and tools during internal transport and during cleaning/sterilization in various dish disinfectors and autoclaves.

The product is a Medical device of risk class I (MDR EU 2017/74)

Material: Stainless steel EN 1.4301 C/C
Grid size: 1.6 mm



INSTRUMENT TRAY WITH LID		
Art.no.	UDI-DI	Measure
705363	735012544030	150x60x45 mm
705362	735012544023	120x60x45 mm
705361	735012544022	60x60x45 mm
705360	735012544021	100x100x60 mm



HARRY HOLMS™
STAINLESS STEEL PRODUCTS SWEDEN

SCISSOR SUPPORT

Support for scissors matching our instrument trays.

Material: Stainless steel EN 1.4301 C/C



SCISSOR SUPPORT

Art.no.	Length	Depth	Wire diameter
705120	250mm	16mm	3mm
705115	340mm	16mm	3mm
705110	450mm	16mm	3,5mm

INSTRUMENT SPRING

Instrument holders suitable for our instrument trays. For fixing and dividing the surgical instruments.

Material: Stainless steel EN 1.4301 C/C



INSTRUMENT SPRING

Art.no.	Length
705105	130mm
705100	190mm



HARRY HOLMS™
STAINLESS STEEL PRODUCTS SWEDEN

ACCESSORIES FOR WIRE BASKETS AND INSTRUMENT TRAYS

The healthcare staff's commitment to the development of our products leads us to useful accessories and details that facilitate daily work. The accessories fit all wire trays and instrument grids in our range.

Material: Stainless steel EN 1.4301 C/C

MARKING PLATE ROUND

Art.no.	Outer diameter	
303000	30mm	Can be marked against a quote



MARKING PLATE RECTANGULAR

2 fastening holes, both short sides rounded

Art.no.	Measure	
306335	63x35mm	Can be marked against a quote



KEY RING

Art.no.	Outer diameter
301905	19,5mm



LABEL HOLDER

Art.no.	Measure
308145	81x45mm



HARRY HOLMS™
STAINLESS STEEL PRODUCTS SWEDEN



Working with people's health is an important assignment. At Mediplast we are in daily contact with dedicated healthcare personnel who do their utmost to deliver better and safer care. These days there is an increasing requirement to work with more patients in less time and with fewer resources.

Mediplast is by your side, always ready to deliver knowledge and innovative solutions with the aim of improving healthcare.

We believe in close relationships. By being present in various care environments on a daily basis, we are able to see your needs and understand your wishes. And through a high level of accumulated expertise and years of experience, we can provide knowledge and advice that really makes a difference, always with the ultimate aim of providing patients with the best possible care.

A close partner to Nordic healthcare

Mediplast is a Swedish supplier of medical technology products, with the Nordic region as our main market. Our product range consists of our own high-quality products, often developed in close collaboration with our healthcare partners, as well as selected products from well-known suppliers

Mediplast has an extensive product range in medical technology, which contains everything from disposable products to advanced equipment of an investment nature. We offer both high quality in-house developed products and are distributors of a number of well-known brands from world-leading manufacturers. Our products cover areas such as thorax, neurology, orthopedics, radiology, cardiology, wound care, ICU, ENT and stoma.