



danumed® GastroTube



New Optimised Product Design

Gastrostomy probe with retention balloon for long-term enteral nutrition

For superior patient comfort and even better ease of use:

New! Innovative rotating ring ENSwivel®, for easy connection and release of ENFit® connections, without twisting the tubes.

New! 90° danuPlate – soft external retention plate with 90° angulation for secure and discreet fixation.

New! danuClip – tube clamp, prevents reflux of nutrition/stomach contents*.

*Available as a separate product, please see REF on back page.



ENSwivel



90° danuPlate



danuClip*

COMFORT


- **Highest patient comfort** due to soft, flexible silicone
- **Rotating ring ENSwivel®** on all ENFit® connectors
- **90° danuPlate** – structured design allows excellent air circulation in the stoma site area
- **danuClip** – practical, simple and quick to use*
- **Extensive selection of sizes** from 12 Fr to 24 Fr for the individual and perfect fit
- Balloon fill valve with Luer connector, with **colour coding** according to Fr size
- **Length marking** in cm along the tube
- **Display of the balloon volume** on the shaft
- **Convenient accessories** with every GastroTube: a pack of **lubricant** for easy insertion and a **danumed® balloon syringe 10ml** to inflate the retention balloon



SAFETY

- Gastric replacement feeding tube made of **high-quality medical grade silicone**
- **Free of Latex, PVC, DEHP and BPA**
- **Flat rounded tip** allows atraumatic insertion
- **Recessed tip** to prevent trauma
- **Inner retention balloon with high balloon stability** for long usage, secure fixation and to avoid leakage
- **ENFit®** connectors according to DIN EN ISO 80369-3 avoid tubing misconnections with other systems
- Available with **ENFit®** and **Y-ENFit®** connectors
- **Radiopaque stripe** along the shaft

Easy handling and secure, discreet fixation of the GastroTube with the 90° danuPlate:

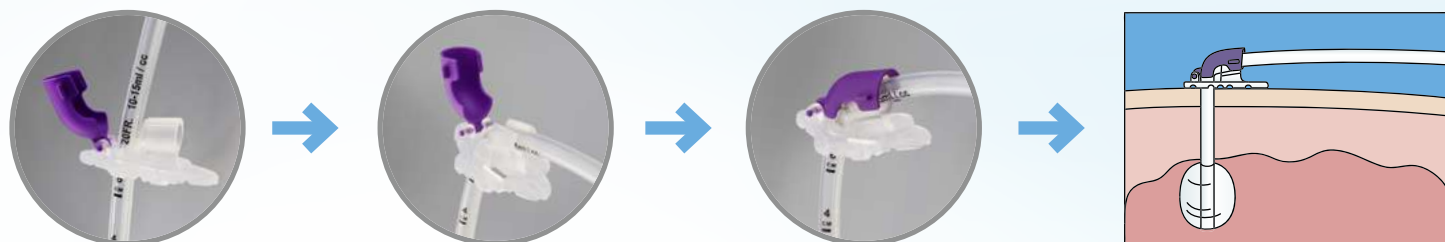

danumed® GastroTubes with ENFit® and Y-ENFit® Supplied in a set with lubricant and danumed® balloon syringe 10ml.



REF	Product description	Size in Fr	Pieces/sales unit	Latex-free	PVC-free	DEHP-free	Bisphenol-A-free
DA 1044 0112	danumed® GastroTube with ENFit® and ENSwivel®	12	1 piece	●	●	●	●
DA 1044 0114	danumed® GastroTube with ENFit® and ENSwivel®	14	1 piece	●	●	●	●
DA 1044 0116	danumed® GastroTube with ENFit® and ENSwivel®	16	1 piece	●	●	●	●
DA 1044 0118	danumed® GastroTube with ENFit® and ENSwivel®	18	1 piece	●	●	●	●
DA 1044 0120	danumed® GastroTube with ENFit® and ENSwivel®	20	1 piece	●	●	●	●
DA 1044 0122	danumed® GastroTube with ENFit® and ENSwivel®	22	1 piece	●	●	●	●
DA 1044 0124	danumed® GastroTube with ENFit® and ENSwivel®	24	1 piece	●	●	●	●
DA 1044 0212	danumed® GastroTube with Y-ENFit® and ENSwivel®	12	1 piece	●	●	●	●
DA 1044 0214	danumed® GastroTube with Y-ENFit® and ENSwivel®	14	1 piece	●	●	●	●
DA 1044 0216	danumed® GastroTube with Y-ENFit® and ENSwivel®	16	1 piece	●	●	●	●
DA 1044 0218	danumed® GastroTube with Y-ENFit® and ENSwivel®	18	1 piece	●	●	●	●
DA 1044 0220	danumed® GastroTube with Y-ENFit® and ENSwivel®	20	1 piece	●	●	●	●
DA 1044 0222	danumed® GastroTube with Y-ENFit® and ENSwivel®	22	1 piece	●	●	●	●
DA 1044 0224	danumed® GastroTube with Y-ENFit® and ENSwivel®	24	1 piece	●	●	●	●
DA 1044 0700	danuClip Tube Clamp	--	15 pcs.	●	●	●	●

Each danumed® GastroTube is individually sterile packed. The danumed® GastroTube is delivered in a sterile peel-pouch and additionally with lubricant and a sterile danumed® balloon syringe 10ml in each single box. 10 boxes in a shipping carton.

Related products to danumed® GastroTubes:



danumed®
Enteral Syringes,
Single Use and
Single-Patient Use



danumed® Gravity Spike,
Gravity Universal and
Gravity Feedbag
Gravity feeding sets

Forhandler:
MEDI PLAST
Medioplast AS
Bjørnstjerne Bjørnsonsgate 110
3044 Drammen
T 32 88 11 00
mediplast@mediplast.no
www.mediplast.com

danumed Medizintechnik GmbH – Individual solutions are our standard.

We are an innovative, German medical technology company with the highest quality standards, many years of expertise in enteral nutrition and clear customer focus. We also meet individual customer requests with flexibility and reliability. Personal, competent, collaborative.

danumed Medizintechnik GmbH • Stadthof 18 • 93059 Regensburg • Germany
Tel + 49 (0) 941/ 890 588-80 • Fax + 49 (0) 941/ 890 588-81
info@danumed.com • www.danumed.com



INNOVATORS IN ENTERAL CARE