




# danumed® StomaMeter

## Stoma Measuring Device with retention balloon

*The StomaMeter is a measuring instrument for determining the stoma length. The length of the stoma tract is measured simply, accurately and reliably to ensure the perfect selection of the optimal gastrostomy button.*

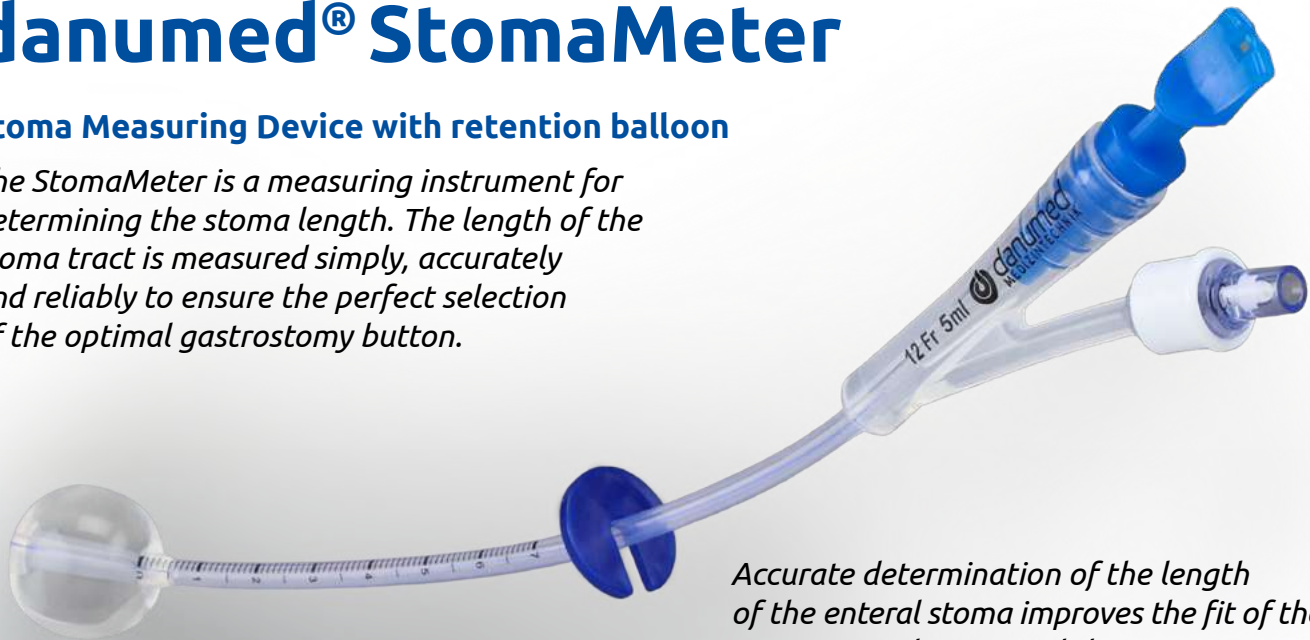

*Accurate determination of the length of the enteral stoma improves the fit of the gastrostomy button and thus patient safety and comfort.*

## COMFORT

- **Soft, flexible material** for best patient comfort
- **Non-slip retention plate** simplifies reading of measured values
- **Imprint of the balloon volume** on the funnel
- **Convenient accessories** in every single package: A pack of **lubricant** for easy insertion and a **danumed® balloon syringe 10ml** to inflate the retention balloon.
- Available in **size 12 Fr**, suitable for each stoma

## SAFETY

- Made of high-quality **medical grade silicone**
- **Latex-free, PVC-free, DEHP-free** and **BPA-free** material
- **Flat rounded tip** to avoid trauma
- **Radiopaque stripe** along the shaft
- **Precisely readable** and wipe-resistant scaling



## danumed® StomaMeter

REF	Size in French	Pieces/sales unit
DA 1045 0112	12	1 piece

The danumed®StomaMeter stoma measuring device is individually sterile packed. The StomaMeter is delivered in a sterile peel-pouch and additionally with lubricant and a sterile danumed® balloon syringe 10ml in each single box. 10 boxes in a shipping carton.

The danumed®StomaMeter is

- **Latex-free • PVC-free • DEHP-free • BPA-free**

Forhandler:  
**MEDI PLAST**  
Medioplast AS  
Bjørnstjerne Bjørnsonsgate 110  
3044 Drammen  
T 32 88 11 00  
mediplast@mediplast.no  
www.mediplast.com