



ORDERING INFORMATION

Adelante® Magnum Set - Includes (1) One Sheath, and (1) One Dilator

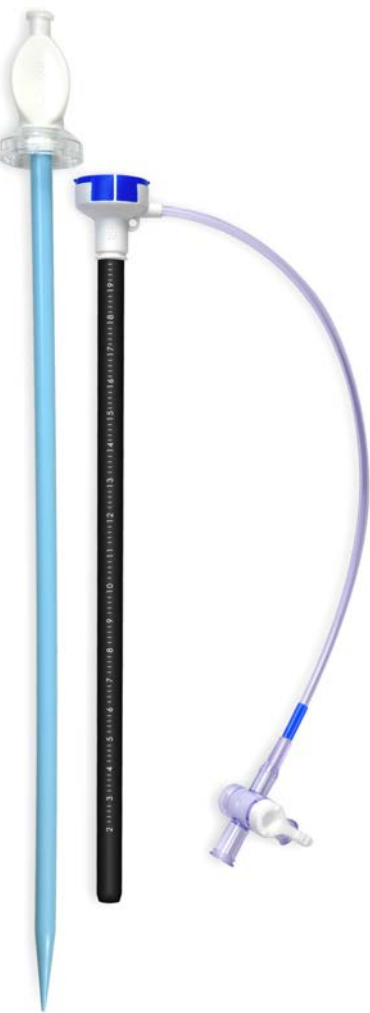

MODEL # SET <i>(One per box)</i>	FRENCH SIZE	SHEATH LENGTH	DILATOR LENGTH
16 F			
ASM016020	● 16 F	20 cm	25 cm
ASM016030	● 16 F	30 cm	35 cm
ASM016035	● 16 F	35 cm	40 cm
18 F			
ASM018020	● 18 F	20 cm	25 cm
ASM018030	● 18 F	30 cm	35 cm

Introducers supplied sterile. For single use only. Available bulk upon request.
Standard lengths are noted above. Other lengths available upon request.
For additional information, e-mail us at sales@oscor.com

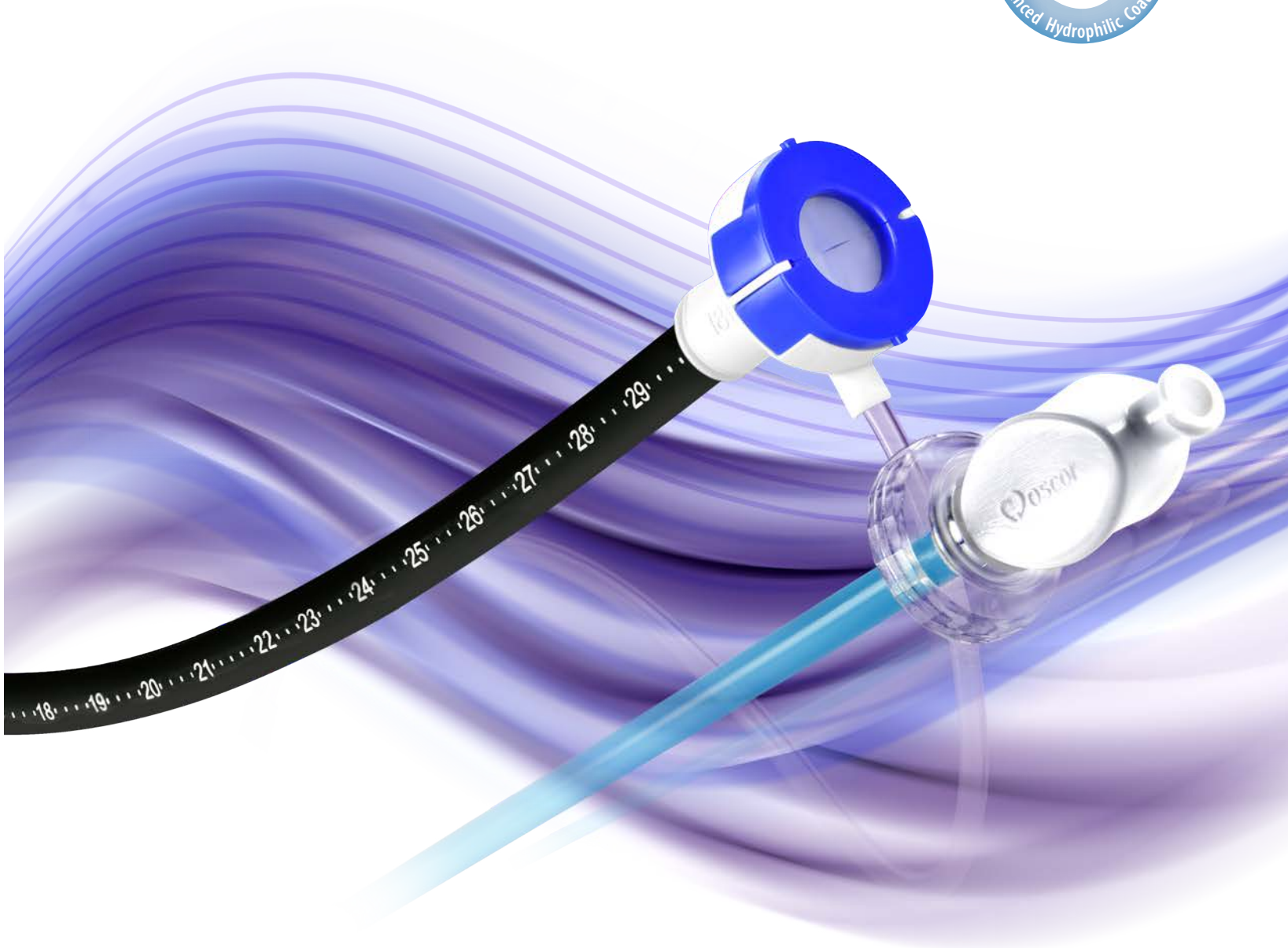

Adelante® **MAGNUM**

HEMOSTATIC VALVE INTRODUCER
FOR LARGE VASCULAR ACCESS

Återförsäljare:
MEDI PLAST
Mediplast AB
Box 1004, 212 10 Malmö
T 020-78 80 35
mediplast.info@mediplast.com
www.mediplast.com

Forhandler:
MEDI PLAST
Mediplast A/S
Marielundvej 46E, 2730 Herlev
T 43 44 40 00
info.dk@mediplast.com
www.mediplast.com

Forhandler:
MEDI PLAST
Mediplast AS
Tollbugata 115, 3041 Drammen
T 32 88 11 00
mediplast@mediplast.no
www.mediplast.com

Jälleenmyyjä:
MEDI PLAST
Mediplast Fenno Oy
PL 153 , 01531 Vantaa
Puh : 09 276 360
info@mediplastfenno.fi
www.mediplast.com
www.mediplastfenno.fi

Återförsäljare:
MEDI PLAST
Mediplast AB
Box 1004, 212 10 Malmö
T 020-78 80 35
mediplast.info@mediplast.com
www.mediplast.com

Forhandler:
MEDI PLAST
Mediplast A/S
Marielundvej 46E, 2730 Herlev
T 43 44 40 00
info.dk@mediplast.com
www.mediplast.com

Forhandler:
MEDI PLAST
Mediplast AS
Tollbugata 115, 3041 Drammen
T 32 88 11 00
mediplast@mediplast.no
www.mediplast.com

Jälleenmyyjä:
MEDI PLAST
Mediplast Fenno Oy
PL 153 , 01531 Vantaa
Puh : 09 276 360
info@mediplastfenno.fi
www.mediplast.com
www.mediplastfenno.fi



P: 727.937.2511
F: 727. 934.9835
E: sales@oscor.com
www.oscor.com

CAUTION: Federal Law restricts this device to sale by or on the order of a physician (U.S.A). Please refer to the Instructions for Use for a complete listing of the indications, contraindications, precautions and warnings, where applicable.



AUG 2014

The **Adelante® MAGNUM** is a high performance hemostatic valve introducer for large size vascular access, specifically designed and optimized for the introduction and placement of endovascular catheters, aortic valves, and stent graft systems. The Adelante® Magnum features a proprietary multilayer hydrophilic coated sheath, having a special designed distal tip portion to accept and tolerate the insertion and manipulation of large and complex devices. The Adelante® Magnum SureSeal™ technology facilitates optimum hemostasis, allowing the insertion of multiple catheters and wires simultaneously with minimal blood loss.

Features:

- Reliable SureSeal™ hemostatic valve technology to minimize blood loss and air embolism even with the insertion of several devices and wires
- Fluid-activated hydrophilic coated sheath, featuring Oscor's proprietary EnSilk™ coating, for exceptional sheath insertion and access
- Multi-layer sheath design provides exceptional kink resistance, while maintaining superb flexibility
- Clear markings on sheath to indicate insertion depth
- Enhanced atraumatic distal sheath tip with radiopaque marker band for accurate positioning under fluoroscopy
- Dilator locks secure into sheath and accepts maximum guidewire up to 0.038"/0.96 mm

