


Safe-Turn Pro & Airway Cradle and Mirror Base



- The unique hand grip on the Cradle give you better control when turning a patient for surgery.
- Adjustable 3 prong Mirror Base with easy to adjust cogs to help achieve a neutral neckline.
- Each Cradle is supplied with both sizes of Prone Face Cushion and a Head Measuring Band.
- View your patient's endotracheal and nasal tubes and their eyes in the adjustable mirror platform.

- Use the tubing clips to keep the tubing in place and stop it tangling.
- Use the head strap a to help you support your patient's head when turning them from supine to the prone position. The head strap should be used as an aid to your usual turning procedure, not as a substitute.
- Although the Safe-Turn is provided with a mirror base, the Cradle can be used with the mirror base to create a lower position for the patient. An additional mirror is provided to be placed under the Cradle.

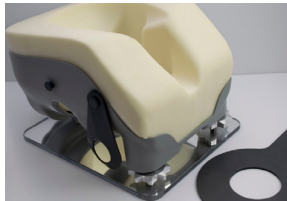

REF ST-PRO-01



REF

Safe-Turn Pro: ST-PRO-01
Box Quantity: 1 pc per box
Box size: 350 x 300 x 300 (mm)

Each box includes;
x1 Safe-Turn Pro Cradle
x1 Safe-Turn Pro Mirror Base with cogs
x2 Safe-Turn Clips
x1 Safe-Turn Strap
x1 Small Mirror (210x130mm)
x1 Prone Face Cushion, Large
x1 Prone Face Cushion, Small
x1 Head Measuring Band
x1 Storage Bag
Plus individual and box barcode labels.

+44 1704 544 944

 **Absolute BPC**

sales@absolutebpc.com

Copyright © 2021. All Rights Reserved.

Safe-Turn Pro & Airway

Cradle and Mirror Base



Safe-Turn Airway Cradle



Safe-Turn Pro Cradle

The Safe-Turn Airway offers all the same features of the Safe-Turn Pro with one difference. The extra cut in the Cradle is designed to accommodate more rigid ET Tubing. The Safe-Turn Airway also allows you to position a patient without the need to disconnect the patient's ET Tubing (when using the Split Prone Face Cushion).

Distributed by:

Återförsäljare:

MEDI PLAST

Mediplast AB

Box 1004, 212 10 Malmö

T 020-78 80 35

mediplast.info@mediplast.com

www.mediplast.com



Easy cleaning for improved infection control. Clean the complete Safe-Turn System with standard hospital wipes.
(Do not Autoclave, Do not submerge in water).

REF ST-AIRWAY-01



REF

Safe-Turn Pro: ST-AIRWAY-01

Box Quantity: 1 pc per box

Box size: 350 x 300 x 300 (mm)

Each box includes;

x1 Safe-Turn Airway Cradle

x1 Safe-Turn Pro Mirror Base with cogs

x2 Safe-Turn Clips

x1 Safe-Turn Strap

x1 Small Mirror (210x130mm)

x1 Prone Face Cushion, Large

x1 Prone Face Cushion, Small

x1 Head Measuring Band

x1 Storage Bag

Plus individual and box barcode labels.

+44 1704 544 944

 **Absolute BPC**

sales@absolutebpc.com

Copyright © 2021. All Rights Reserved.