



## SUCTION CATHETER


Straight, coudé (anged tip)  
as well as with vacuum control

Medioplast has one of Scandinavia's biggest assortments of high quality suction catheters, manufactured for all the different needs in care. Dependent on the area of use, we not only offer various lengths and sizes but also various designs, materials and packaging.


## DESIGN

All suction catheters have a central, soft hole at the tip **without any sharp edges**. There are also two opposite relief holes, about one cm from the end to secure an atraumatic suction.


The coudé catheter has an **anged tip** that makes it possible to control the direction of the catheter in the bronchial tree. A little lump on the connector indicates the direction for which the tip is bent.

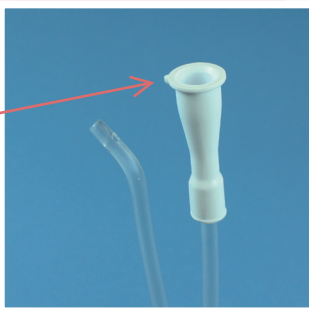


Suction catheters with built-in **vacuum control** are available. These will enable you to easily control the suction strength with your fingertip through intermittent suction technique.

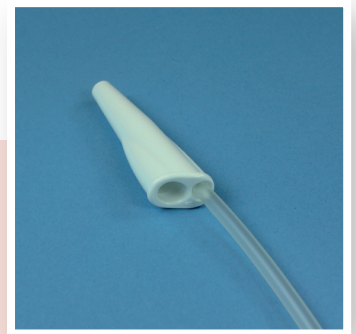


All suction catheters have a colour coded connector according to the standard for suction catheters **EN ISO 8836:2019** in order to facilitate identification of the correct size.

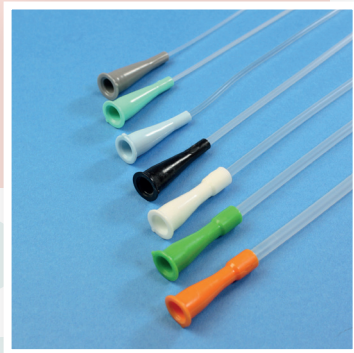


Straight catheters are available in **3 lengths**, incl. connector (4 cm), 25 cm, 40 cm, 53 cm.

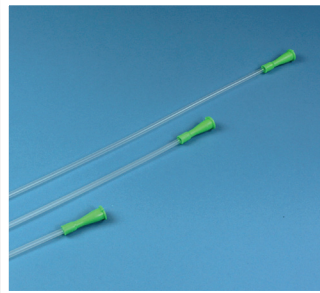

Suction catheters with vacuum control are available in **two lengths**, incl. connector (6 cm): 35 cm suction catheter with centimeter graduation as well as 53 cm without graduation.

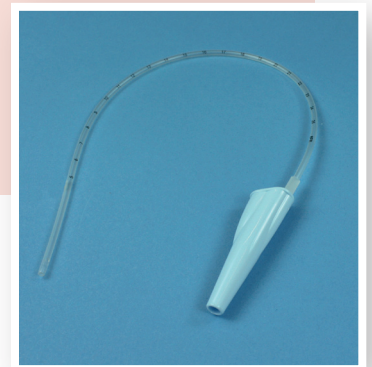

Coudé suction catheter (anged tip) are available in 53 cm length, incl. connector (4 cm).

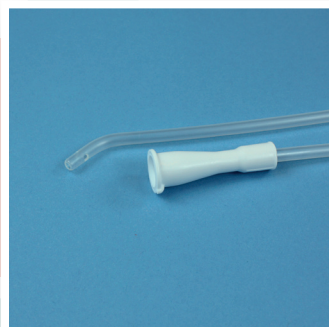

All suction catheters are single-use, delivered sterile and with a **shelf life of 5 years**. Sterilisation method, ethylene oxide.



## MATERIALS


Mediplast offers suction catheters in two different materials, **PVC without phthalates or TPE** which is a PVC-free material. The softening agent (plasticizer) used for the PVC is a non-phthalate, called ATBC.

The material in a suction catheter is rarely recycled, and in most cases it is incinerated. Because of environmental concerns, today there is an effort to not use PVC, as it consists of chlorinated plastic that forms hydrochloric acid (HCl) during combustion. PVC requires specialized facilities to handle the hydrochloric acid that is formed, and since there are PVC-free alternatives, they are often recommended in stead.

For the user, there could be a small, insignificant difference in the two materials. The PVC-free alternative is a little more stable, whereas the suction catheter of PVC is a little softer and more flexible. The difference is hardly detectable but when the area of use varies, the choice of material may also be of significance.

## PACKAGING


The product packaging is a so-called **peel pack** made of paper and plastic. It has a hole at the top to hang it.

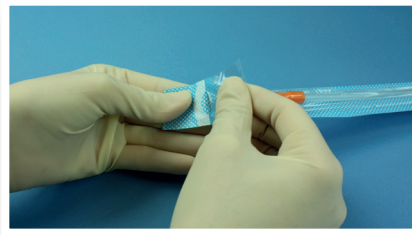


The packaging also has a blue dotted interior, which is a coating that makes it easy to separate the paper from the plastic **without damaging it**. A safe and aseptic opening is secured.

The products are packaged according to the three-layer principle. 100 pc of suction catheters in one carton ward box, 6 ward boxes with 600 suction catheters in one transport box. The smallest unit that can be ordered is 100 pc = MBE.

To simplify the removal from the ward box, Mediplast is now able to offer **individually packaged** suction catheters. The product is easier to remove, while at the same time you can be certain that the catheters are kept straight and flat inside the box.

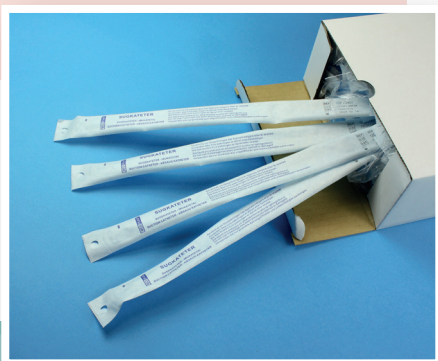

The standard packaging is different as the suction catheters are placed in **perforated sheets** in the carton ward box. This still enables the removal of one suction catheter at a time but it simplifies the distribution in those circumstances where you want to place the catheters in a holder.

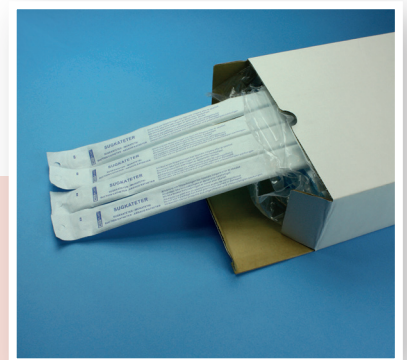

**SUCTION CATHETER STRAIGHT, 25 CM**

PVC ATBC single pack	PVC-free single pack in perforated sheets	PVC ATBC single pack in perforated sheets	Size (Ch)	Colour	Ward/Trsp box
	60880525	60130525-1	05	grey	100/600
608806251	60880625	60130625-1	06	light green	100/600
608808251	60880825	60130825-1	08	light blue	100/600
608810251	60881025	60131025-1	10	black	100/600
608812251	60881225	60131225-1	12	white	100/600
608814251	60881425	60131425-1	14	green	100/600
	60881625		16	orange	100/600


**SUCTION CATHETER STRAIGHT, 40 CM**

PVC ATBC single pack	PVC-free single pack in perforated sheets	PVC ATBC single pack in perforated sheets	Size (Ch)	Colour	Ward/Trsp box
	60880540		05	grey	100/600
608806401	60880640	60130640-1	06	light green	100/600
608808401	60880840	60130840-1	08	light blue	100/600
608810401	60881040	60131040-1	10	black	100/600
608812401	60881240	60131240-1	12	white	100/600
608814401	60881440	60131440-1	14	green	100/600
608816401	60881640	60131640-1	16	orange	100/600

**SUCTION CATHETER STRAIGHT, 53 CM**


PVC ATBC single pack	PVC-free single pack in perforated sheets	PVC ATBC single pack in perforated sheets	Size (Ch)	Colour	Ward/Trsp box
608806531	60880653	60130653-1	06	light green	100/600
608808531	60880853	60130853-1	08	light blue	100/600
608810531	60881053	60131053-1	10	black	100/600
608812531	60881253	60131253-1	12	white	100/600
608814531	60881453	60131453-1	14	green	100/600
608816531	60881653	60131653-1	16	orange	100/600
		60131853-1	18	red	100/600

Ward box = min qty / order

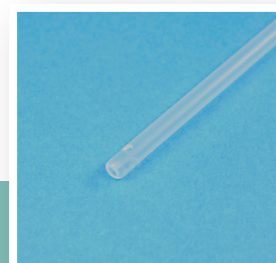

## SUCTION CATHETER, COUDÉ, 53 CM

PVC ATBC single pack	PVC-free single pack in perforated sheets	PVC ATBC single pack in perforated sheets	Size (Ch)	Colour	Ward/ Trsp box
608708531	60870853	60170853-1	08	light blue	100/600
608710531	60871053	60171053-1	10	black	100/600
608712531	60871253	60171253-1	12	white	100/600
608714531	60871453	60171453-1	14	green	100/600
608716531	60871653	60171653-1	16	orange	100/600

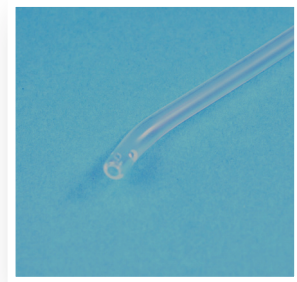

## SUCTION CATHETER WITH VACUUM CONTROL, CM-GRADUATION, 35 CM

PVC ATBC - single pack in perforated sheets	Size (Ch)	Colour	Ward/ Trsp box
60151804-1	04	purple	100/600
60151805-1	05	grey	100/600
60151806-1	06	light green	100/600
60151808-1	08	light blue	100/600
60151810-1	10	black	100/600

## SUCTION CATHETER WITH VACUUM CONTROL, 53 CM

PVC ATBC - single pack in perforated sheets	Size (Ch)	Colour	Ward/ Trsp box
60150653-1	06	light green	100/600
60150853-1	08	light blue	100/600
60151053-1	10	black	100/600
60151253-1	12	white	100/600
60151453-1	14	green	100/600
60151653-1	16	orange	100/600
60151853-1	18	red	100/600

Ward box = min qty / order



Working with people's health is an important assignment. At Mediplast we are in daily contact with dedicated healthcare personnel who do their utmost to deliver better and safer care. These days there is an increasing requirement to work with more patients in less time and with fewer resources.

Mediplast is by your side, always ready to deliver knowledge and innovative solutions with the aim of improving healthcare.

We believe in close relationships. By being present in various care environments on a daily basis, we are able to see your needs and understand your wishes. And through a high level of accumulated expertise and years of experience, we can provide knowledge and advice that really makes a difference, always with the ultimate aim of providing patients with the best possible care.

A close partner to Nordic healthcare.

Mediplast is a Swedish supplier of medical technology products. Our range consists of high-quality products under our own brand as well as selected products from well-known suppliers.

#### **Disposal**

Peelpack: Divided into two parts - Paper → paper disposal

Plastic (PP/PE) → plastic disposal

Inner pouch in wardbox: LDPE → plastic disposal

Wardbox, white: → paper disposal

Transportbox, brown: → paper disposal