

TWIST

 **oscor**® | *for life!*®



Destino™ **TWIST**

DEFLECTABLE STEERABLE GUIDING SHEATH

C: 9 mm

 **oscor**®

6.5 F
UL: 45 cm
OL: 60 cm



The *Destino™ TWIST* is a steerable guiding sheath that offers a deflectable tip with an ergonomic handle. The Destino™ Twist is the ultimate tool for gaining access to the most difficult to reach sites, while maintaining hemostasis with Oscor SureSeal™ technology. The Destino™ Twist is **YOUR ONE SHEATH SOLUTION** intended to facilitate the intracardiac, renal, and peripheral placements of diagnostic and therapeutic devices.

Clear view of the distal marker band under fluoroscopy

Smooth tip to dilator transition

Hydrophobic coated sheath

Braided Flexsteer™ shaft technology

Rotating collar

Tip position indicator

Reliable SureSeal™ hemostatic valve

French size and guidewire indicator

6.5F .035"

C: 9 mm

6.5 F
UL: 45 cm
OL: 60 cm

0


0
C: 9 mm

INTERVENTIONAL USE



SHAFT MODELS


Tichà Damascelli (TD) Curve

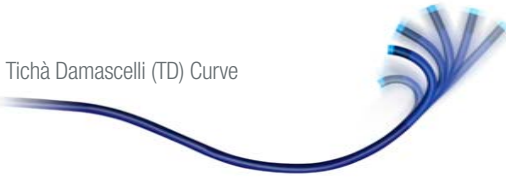

Standard



SHAFT DIAMETERS

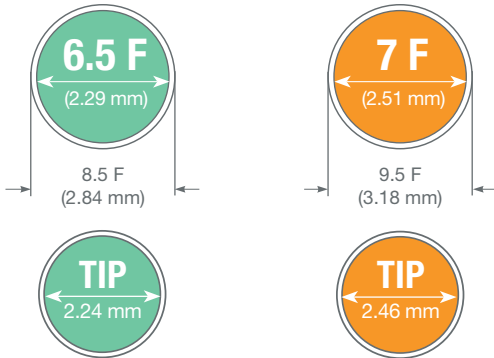

TIP DEFLECTION



TRANSSEPTAL USE




SHAFT MODELS


Standard



SHAFT DIAMETERS

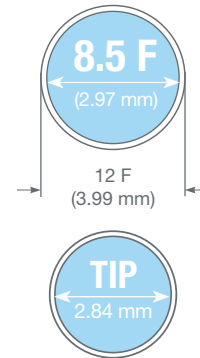

TIP DEFLECTION

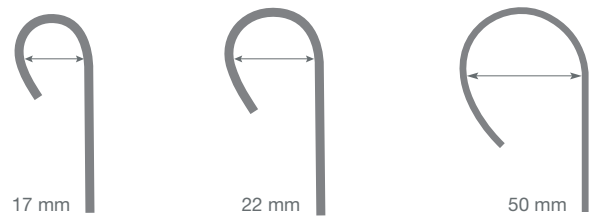


Features:

- Precise tip deflection adjustment for fast access
- Tip position indicator showing neutral or deflected curve
- Superior curve retention for reliable catheter placement
- Enhanced radiopaque tip visibility enabling accurate tip positioning
- Proprietary hydrophobic coated sheath technology, enabling smooth access
- Reliable Flexsteer™ reinforced braid technology provides exceptional control, kink resistance, while maintaining optimal flexibility
- Ergonomic handle with clear imprint of the curve deflection, french size and sheath length
- Reliable SureSeal™ hemostatic valve technology to minimize air embolism and blood loss even under extreme conditions



ORDERING INFORMATION

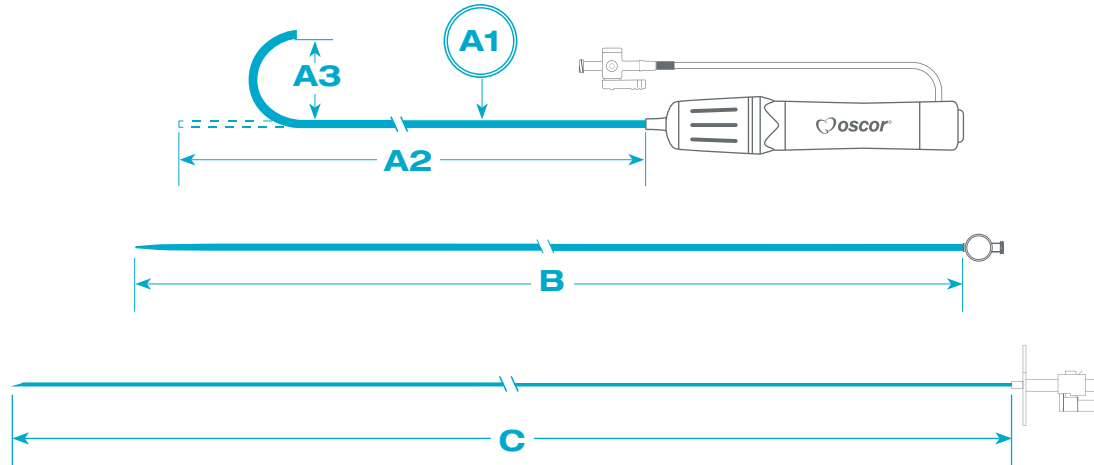

INTERVENTIONAL USE - ORDERING INFORMATION

Model # Set <i>(One per box)</i>	A1 Sheath Inner Diameter	A2 Usable Sheath Length	B Usable Dilator Length	A3 Curve Bend
STANDARD				
DST0654509	● 6.5 F	45 cm	63 cm	9 mm
DST0655509	● 6.5 F	55 cm	73 cm	9 mm
DST0659009	● 6.5 F	90 cm	108 cm	9 mm
DST0704509	● 7 F	45 cm	63 cm	9 mm
DST0705509	● 7 F	55 cm	73 cm	9 mm
DST0705517	● 7 F	55 cm	73 cm	17 mm
DST0709009	● 7 F	90 cm	108 cm	9 mm
TD CURVE				
DST0654509TD	● 6.5 F	45 cm	63 cm	9 mm
DST0655509TD	● 6.5 F	55 cm	73 cm	9 mm
DST0704509TD	● 7 F	45 cm	63 cm	9 mm
DST0705509TD	● 7 F	55 cm	73 cm	9 mm

TRANSSEPTAL USE - ORDERING INFORMATION

Model # Set <i>(One per box)</i>	A1 Sheath Inner Diameter	A2 Usable Sheath Length	B Usable Dilator Length	A3 Curve Bend	C Compatible Transseptal Needle
STANDARD					
DST0856717	● 8.5 F	67 cm	85 cm	17 mm	89 cm
DST0857317	● 8.5 F	73 cm	91 cm	17 mm	98 cm
DST0856722	● 8.5 F	67 cm	85 cm	22 mm	89 cm
DST0857322	● 8.5 F	73 cm	91 cm	22 mm	98 cm
DST0856750	● 8.5 F	67 cm	85 cm	50 mm	89 cm
DST0857350	● 8.5 F	73 cm	91 cm	50 mm	98 cm

Guiding sheaths are supplied sterile. For single use only. Available bulk upon request.
Standard lengths are noted above. Other lengths available upon request.



P: 727.937.2511
F: 727.934.9835
E: sales@oscor.com
www.oscor.com

CAUTION: Federal Law restricts this device to sale by or on the order of a physician (U.S.A.). Please refer to the Instructions for Use for a complete listing of the indications, contraindications, precautions and warnings, where applicable.



Återförsäljare:

MEDI PLAST

Mediplast AB
Box 1004, 212 10 Malmö
T 020-78 80 35
mediplast.info@mediplast.com
www.mediplast.com

Forhandler:

MEDI PLAST

Mediplast A/S
Marielundvej 46E, 2730 Herlev
T 43 44 40 00
info.dk@mediplast.com
www.mediplast.com

Forhandler:

MEDI PLAST

Mediplast AS
Tollbugata 115, 3041 Drammen
T 32 88 11 00
mediplast@mediplast.no
www.mediplast.com

Jälleenmyyjä:

MEDI PLAST

Mediplast Fenno Oy
PL 153 , 01531 Vantaa
Puh : 09 276 360
info@mediplastfenno.fi
www.mediplast.com
www.mediplastfenno.fi